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H_LAG3 Reporter Jurkat Cell Line

H_LAG3 Reporter Jurkat 细胞

For research use only!

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Reagent	Specification	Manufacturer/Catalogue No.
Puromycin	25 mg	Genomeditech/GM-040401-1
Blasticidin	10 mg	Genomeditech/GM-040404-1
G418 Sulfate	1 g	Genomeditech/GM-040402-1
金黄色葡萄球菌肠毒素 (SEE 肠毒素)	/	Genomeditech/GM-H23036
Fetal Bovine Serum	500 mL	ExCell/FSP500
RPMI 1640	500 mL	Gibco/C11875500BT
96 Well Clear V-Bottom Tissue Culture	96-Well	Corning/3894
96 Well round Well culture plate	96-Well	NEST/701001
96 Well White Flat Bottom Polystyrene Not Treated Microplate	96-Well	Corning/3912
Switch-On Reagent(1000X)	1 mL	Genomeditech/GM-041519
Anti-H_LAG3 hIgG4 Antibody(Mavezelimab)	/	Genomeditech/GM-28753AB
Anti-Lag3 hIgG4 Antibody(Relatlimab)	/	Genomeditech/GM-51650AB
Anti-H_CD4 hIgG1 Antibody(Tregalizumab)	/	Genomeditech/GM-28752AB
Switch-On	/	Genomeditech/GM-041519
GMOne-Step Luciferase Reporter Gene Assay Kit	1000T	Genomeditech/GM-040503C

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Equipment	Manufacturer/Catalogue No.
细胞计数仪	ThermoFisher Scientific/Countess 3
酶标仪	Moleculardevices/SpectraMax L

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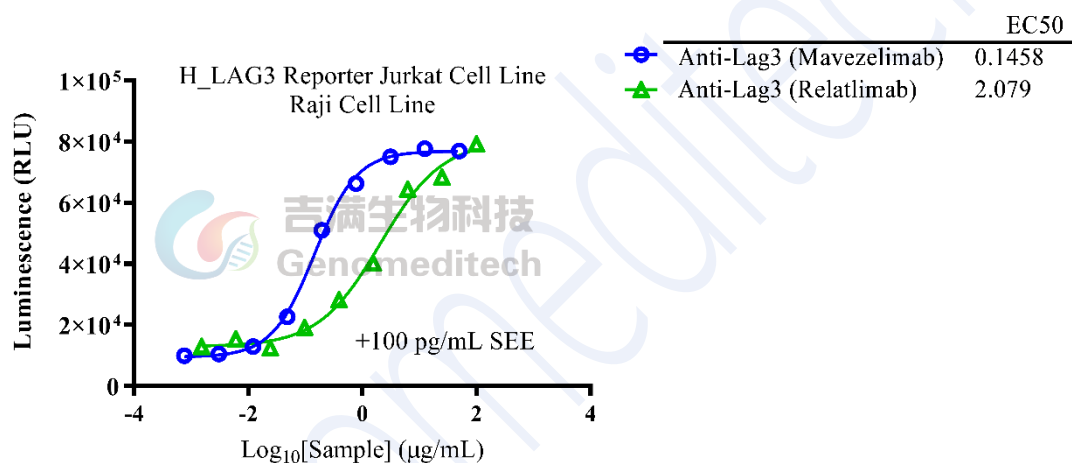
Fig 1. Response to Mavezelimab and Relatlimab. The H_LAG3 Reporter Jurkat Cell Line (Cat. GM-C20096) was pre-treated with Switch-On Reagent(Cat. GM-041519) for 48 hours prior to the assay. Serial dilutions of Mavezelimab (Cat. GM-28753AB) and Relatlimab (Cat. GM-51650AB) were incubated with 1E5 cells/well of the cells in a 96-well plate for 30 minutes. Simultaneously, 10 pg/well of SEE was added to 2E4 cells/well of Raji cells in another 96-well plate and incubated for 30 minutes. The two mixtures were then combined and incubated for an additional 16 hours. Firefly luciferase activity was measured using the GMOne-Step Luciferase Reporter Gene Assay Kit (Cat. GM-040503). The results indicated a maximum blocking fold of approximately [7.8] and [6.2]. Data are shown by drug mass concentration.

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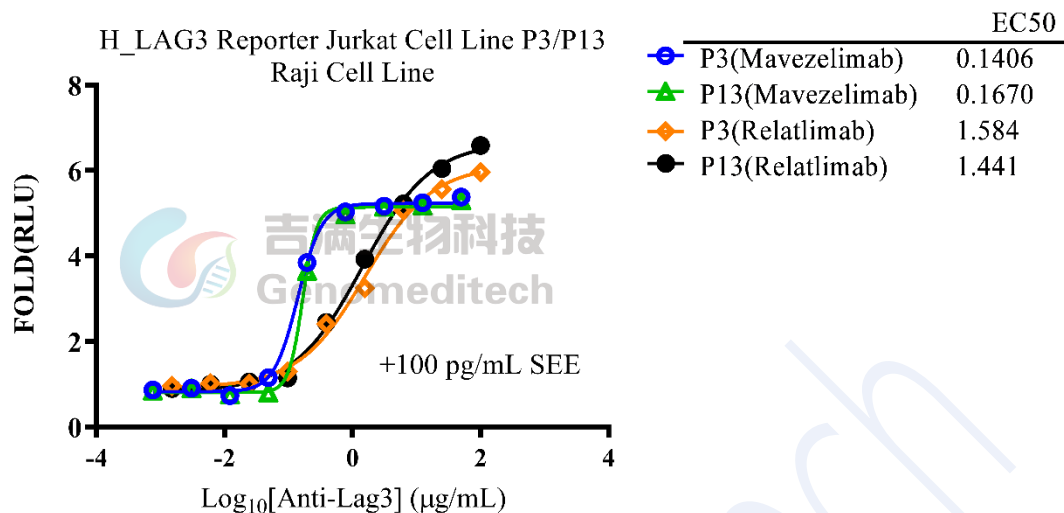
Fig 2. The passage stability of response to Mavezelimab and Relatlimab. The H_LAG3 Reporter Jurkat Cell Line (Cat. GM-C20096) was pre-treated with Switch-On Reagent(Cat. GM-041519) for 48 hours prior to the assay. Serial dilutions of Mavezelimab (Cat. GM-28753AB) and Relatlimab (Cat. GM-51650AB) were incubated with 1E5 cells/well of the passage 3 and 13 of the cells in a 96-well plate for 30 minutes. Simultaneously, 10 pg/well of SEE was added to 2E4 cells/well of Raji cells in another 96-well plate and incubated for 30 minutes. The two mixtures were then combined and incubated for an additional 16 hours. Firefly luciferase activity was measured using the GMOne-Step Luciferase Reporter Gene Assay Kit (Cat. GM-040503). Data are shown by drug mass concentration.

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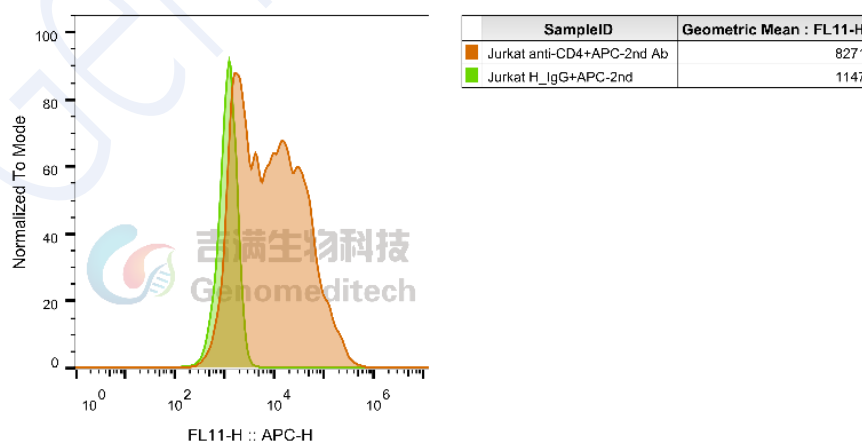
Fig 3. H_LAG3 Reporter Jurkat Cell Line (Cat. GM-C20096) was determined by flow cytometry using Anti-H_CD4 hIgG1 Antibody(Tregalizumab) (Cat. GM-28752AB).

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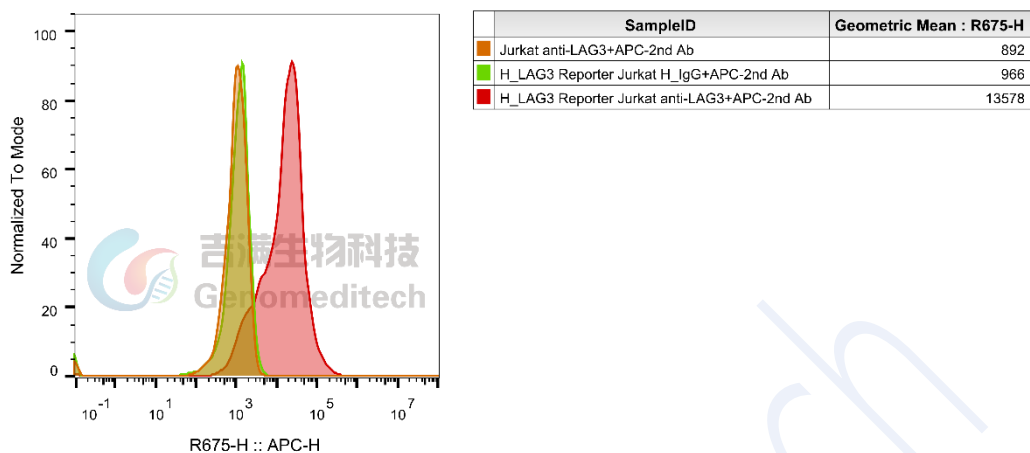
Fig 4. H_LAG3 Reporter Jurkat Cell Line (Cat. GM-C20096) was determined by flow cytometry using Anti-H_LAG3 hIgG4 Antibody(Mavezelimab) (Cat. GM-28753AB) (Cells were induced with Switch-On for 48 h).

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Lag3	
Cynomolgus_LAG3 CHO-K1 Cell Line	H_LAG3 CHO-K1 Cell Line
H_LAG3 HEK-293 Cell Line	Rhesus_LAG3 CHO-K1 Cell Line
Anti-H_LAG3 hIgG4 Antibody(Mavezelimab)	Anti-Lag3 hIgG4 Antibody(IMP-7)
Anti-Lag3 hIgG4 Antibody(Relatlimab)	

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